

PRONE STANDER

PART # 4617A

Box size: 34½" H. x 25" W. x 44½" D. Shipping Weight: 75 lbs.

PATTERSON
MEDICAL

SAMMONS
PRESTON

ADJUSTING TO FIT USER - INITIAL SETUP

1. Remove unit from box. (Fig.1)
2. Turn crank clockwise (Fig.2) to bring stander to near vertical position. (Fig.3)



Fig.1

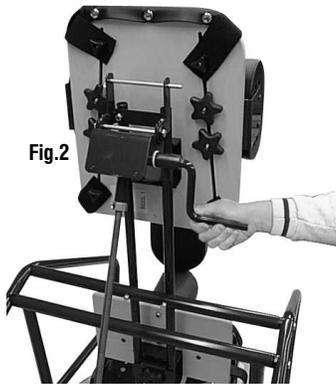


Fig.2



Fig.3

ADJUSTING BODY BOARD, TRUNK SUPPORTS, BODY STRAPS

1. To adjust **Body Board**, securely grasp, and loosen knob. (Fig.4) Adjust vertically to desired position and tighten knob.
2. To adjust **Trunk Supports**, loosen knob, adjust horizontally to desired position and tighten knob. (Fig.5)
3. To adjust **Body Straps**, slide buckles to desired position on straps. **NOTE: The straps must be crossed over users shoulders and back.** (Fig.6)

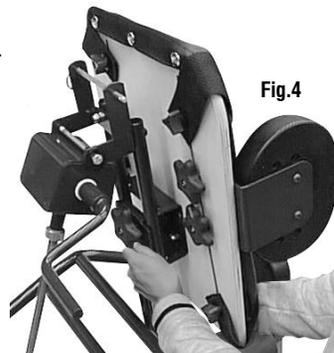


Fig.4

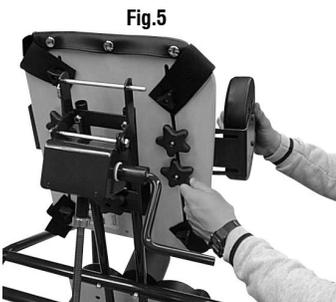


Fig.5



Fig.6

ADJUSTING KNEE BOARD, POSTERIOR SUPPORT

1. To adjust **Knee Board**, securely grasp, and loosen knob. (Fig.7) Adjust vertically to desired position and tighten knob.
2. To adjust **Posterior Support** vertically, securely grasp, and loosen knob. (Fig.8) Adjust to desired position, tighten knob.
3. To adjust **Posterior Support** horizontally, securely grasp, and loosen knob. (Fig.9) Adjust to desired position, tighten knob.

Fig.7



Fig.8



Fig.9

ADJUSTING FOOT BOARD, FOOT HOLDERS, FOOT STRAPS

1. To adjust **Foot Board**, securely grasp, and loosen knob. (Fig.10) Adjust vertically to desired position and tighten knob. **NOTE: Knee Board may be removed, allowing Foot Board to accommodate a shorter user.** (Fig.11)
2. To adjust **Foot Holders**, loosen knobs. Adjust horizontally to desired position and tighten knobs. (Fig.12)

Fig.10

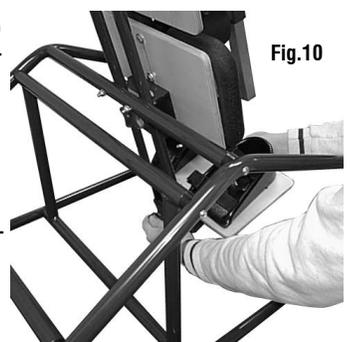


Fig.11

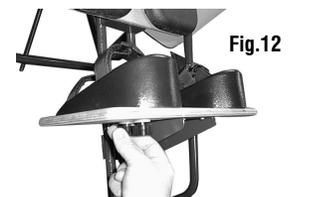


Fig.12

ATTACHING TRAY

1. Remove the ring pin (Fig.14) from the bracket above the crank. (Fig.15)
2. Place the hook bracket over lower tube on front of unit. (Fig.16)
3. Rotate tray assembly into user position, position bracket under tray onto bracket above crank (Fig.17) Connect tray to unit with ring pin.

Fig.14



Fig.15



Fig.16

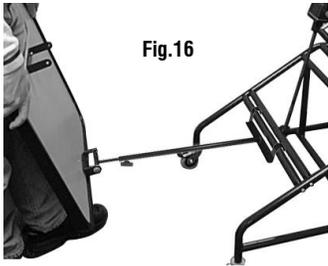


Fig.17




ADJUSTING TRAY ANGLE, STANDER ANGLE

1. To adjust front elevation of **Tray**, loosen knob and slide assembly to desired position. (Fig.18) Tighten knob.
2. To adjust the angle of the stander, turn the crank in either direction until the stander is in the desired position. (Fig.19)

Fig.18



Fig.19

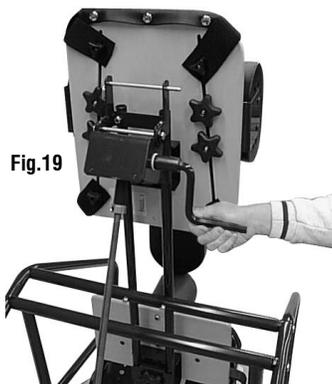

LOCKING WHEELS

1. To lock wheel, push down on black tab. (Fig.20) This snaps tab against wheel, into locked position. (Fig.21)

Fig.20



Fig.21



PRECAUTIONS

- ◆ DO NOT LEAVE CHILD UNATTENDED IN STANDER
- ◆ CONSULT A HEALTHCARE PROFESSIONAL FOR APPROPRIATE CRITERIA FOR INDIVIDUAL USAGE